




Friendly To The Environment - Being green isn't just about the materials, it's also about manufacturing processes, air quality and conservation. Our furniture meets or exceeds stringent manufacturing criteria by controlling air quality and emission levels. Our state-of-the-art computerized UV flat line spray system does not release any volatile organic compounds (VOCs). Our energy efficient manufacturing processes and waste and recycle management is achieved by using the latest wood working machine technology. Using the

newest machinery increases productivity and capacity with less energy. We are certified by the SmartWood Program of the Rainforest Alliance to FSC Chain-of-Custody standards, offer certified and sustainable materials, no added urea-formaldehyde (NAUF), California Air Resources Board (CARB) compliant and we have a LEEP AP (Accredited Professional) on staff ready to assist with LEED® registered projects.

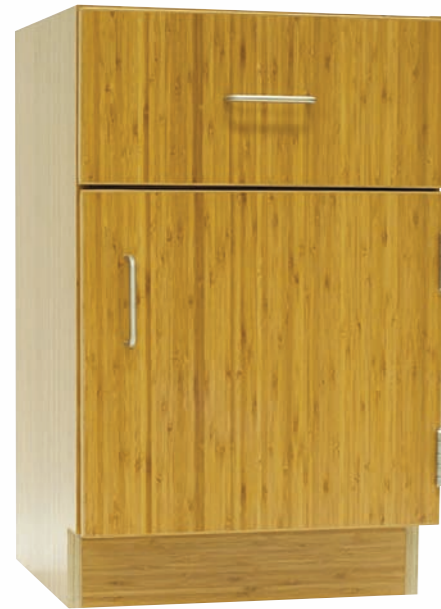

Mott Wood Casework

- Provide FSC certified wood that is supplied by sustainable foresters, harvesters and manufactures which contributes to LEED credit MR7
- Choice of rapidly renewable materials including bamboo which contribute to LEED credit MR6
- Offer NAUF panel products which directly impacts EQ4.4 and are California Air Resources Board (CARB) compliant
- Provide particleboard or fiberboard that has pre-consumer recycled content and recovered material recycled content

For more information on Mott's sustainable solutions; simply contact us with your needs.

QUALITY IN PRODUCT • PRECISION IN PROCESS • STRENGTH IN PEOPLE

WOOD LABORATORY

WOOD LABORATORY FURNITURE



QUALITY BY DESIGN

Mott Manufacturing provides complete solutions for all your laboratory needs. We offer a full line of products featuring quality steel laboratory furniture, premium grade laboratory wood casework, flexible and mobile furniture systems, high performance and custom fume hoods, custom stainless steel work surfaces and sinks. We have been serving the industrial, pharmaceutical, education, health care and government markets since the early 1960's. Since each laboratory's needs are unique, we can design and engineer products to meet your exact requirements.

Customer Satisfaction is Mott's First Priority

Delivering industry-leading customer service is more than just words to the Mott Team. We know that satisfied customers drive our success. We therefore provide our customers with the following key values:

Consistent Product Quality

Sophisticated high-tech machinery, computerized engineering processes, quality materials and skilled craftsmanship combine to assure quality.

Just-in-Time Delivery Capability

Our six week manufacturing lead time, along with our high on-time and complete shipment performance, provide exceptional reliability to our customers.

Adaptable Product Line

Our extensive modular product line permits our customers to readily configure the most productive solution to meet their specific requirements. Our flexible manufacturing capabilities allow convenient modification of standard components to fit the unique requirements of each laboratory.

Competitive Value

Our several decades of application experience allows us to provide productive furniture solutions to our customers. This experience, along with our high product quality, custom manufacturing capability, consistently strong service levels, and competitive pricing provide high value to our customers.

Extensive Dealer Network

Our dealers are the best. We build strong relationships with our dealers to provide, in partnership, industry-leading service levels. Our dealers' on site leadership and years of experience combined with our design and manufacturing skills deliver the very best in laboratory furniture solutions.



Corporate Headquarters
452 Hardy Rd, Brantford ON, Canada N3T 5L8
T. 519.752.7825 • F. 519.752.2895

Wood Casework Division - 562 Industrial Park Rd,
P.O. Box 156, Maxwelton, WV, USA 24957
T. 304.497.2115

www.mott.ca



AUGUST 2011





Quality Craftsmanship....

Traditional quality craftsmanship brings a rich and lasting presence to your environment. Sophisticated high-tech machinery, computerized engineering processes, quality materials and skilled craftsmanship combine to assure quality. With Mott, you are investing in cabinetry engineered with structural integrity, durability and aesthetics proven to withstand the rigors of everyday use.

Our commitment to quality is exceptional. A cabinet is only as good as the craftsmanship that goes into it. We precision manufacture all our components to exceptionally tight tolerances. Computer-cutting capabilities and exceptional finishing ensure that the details of Mott cabinets stands out under the most rigorous scrutiny. All furniture components have been tested to the SEFA 8 Recommended Practices and are built to AWI standards to ensure many years of intensive use in a demanding environment.

Everybody's needs are unique; Mott can design and manufacture to meet your exact requirements. We have the expertise to provide solutions to meet your demands, simply contact us with your needs.



• **Vertical Grain Match** - All Mott cabinets feature vertical grain drawer and door fronts. Optional combination matched fronts an available option.




• **Drawer Body Options** - Drawers are available in Baltic Birch, NAUF 1/2" plywood panels, powder coated cold rolled steel, powder coated galvanized steel and stainless steel.



• **Dovetail Drawers** - Full height drawer sides are constructed with 1/2" thick 9-ply Baltic Birch with durable dovetail joinery at all four corners. Hardboard, melamine or veneered drawer bottoms are available options.



• **Full Frame Support Platform** - Optional full frame construction provides maximum strength and rigidity to allow for cabinets to be easily handled during installation or renovation.

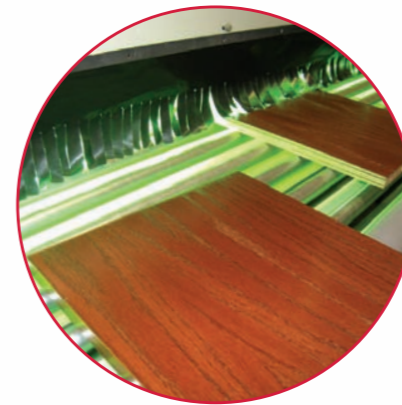

• **UV Flat Line Finishing** - Our flat line finish system offers great chemical resistivity, emits zero volatile organic compounds (VOCs) using solvent-free UV coatings.

- **Cabinet Construction** - 3/4" thick cabinet case panels are joined together with concealed hardwood dowels and NAUF glue, providing structural integrity and strength.
- **Drawer and Door Fronts** - 3/4" particleboard core banded on all sides with edge banding that matches cabinet veneer. MDF and combo core are available options.
- **End Panels** - 3/4" veneer core with hardwood edge banding.
- **Sturdy Bottom** - 3/4" veneer core with a 1/8" hardwood edge band set flush and attached to the cabinet ends with dowels.
- **Drawer Suspension** - Several different drawer suspension systems available including 100# nylon tire self closing, 100# full extension, 150# full extension and 200# full extension for heavy-duty requirements. Full extension provides full access to all corners of the drawer.
- **Adjustable Shelves** - Fully adjustable shelving supported by seismic clips for storage flexibility. All shelves are of the same species as exterior and have 1/8" hardwood edgeband on front edge. Shelves are adjustable in 1-1/4" increments for maximum flexibility.
- **Quality Hardware** - Five knuckle, three knuckle semi-concealed and self-closing European style hinges and a wide selection of handles available to meet specific aesthetic and functional requirements.
- **Exterior Options** - We can provide a wide range of wood veneers, in various cuts, matching and sequencing options. We offer a variety of standard stains and custom staining to match existing casework is available.
- **Full Range of Cabinetry and Components** - For all of your needs, our base, wall and floor cabinets are available in numerous configurations, in a wide variety of heights and widths with countless options.
- **Mobility** - All of our base cabinets and tables are available with heavy-duty locking casters providing convenient storage and mobility allowing users to quickly and easily move to another workstation as tasks change.
- **Student Workstations** - A wide selection of Student Table Assemblies designed to provide flexibility and durability in classroom laboratory environments are available. Choose from island and wall assemblies available in numerous configurations and sizes with countless options. Complete the assemblies with sinks, electrical, plumbing fixtures and accessories. Consult with a Mott dealer representative to discuss your specific requirements.

Full Flush Overlay



- Features square-edged flush overlay door and drawer fronts with 1/16" vertical reveals at the cabinet ends and 1/8" vertical and horizontal reveals between doors/drawers.

Reveal Overlay



- Features square-edged door and drawer fronts, with 1/4" vertical reveals at the cabinet ends and 1/8" vertical and horizontal reveals between doors and drawers.

Lipped Radius Overlay



- Features radiused edged door and drawer fronts machined to be inset 3/8" from the face of the cabinet, with 1/2" vertical reveals at the cabinet ends and 1/8" vertical and horizontal reveals between doors/drawers. Cabinet fronts have a combination grain direction.